

## The Pursuit of Safety throughout Infusion

- **Dynamic Pressure Sensors (DPS)**  
Real-time display of pressure changes in the pipeline, caregivers can predict the risk through obstruction early warning to ensure the safety of drug injection.
- **Automatic pressure release (Anti-Bolus)**  
When the pipeline is blocked and alarms, the pressure in the pipeline is automatically released to avoid unexpected bolus to the patient.



## Inter-connection Empowers Imagination

Support wired / WiFi connectivity to meet stable transmission for clinical information.  
Access to MX8900 Infusion Workstation System.  
Access to HIS/CIS/EMR through HL7 protocol.

## M800 Syringe Pump




**Shenzhen Comen Medical Instruments Co., Ltd.**

Add: Floor 10, Floor 11 and Section C of Floor 12 of Building 1A & Floor 1 to Floor 5 of Building 2, FIYTA Timepiece Building, Nanhuan Avenue, Matian Sub-district, Guangming District, Shenzhen, Guangdong, 518106, P.R. China  
Tel: +86-72643 1236 Fax: +86-755-2643 1232 Website: en.comen.com E-mail: info@szcomen.com

P/N: EN-M800-4P-20211026-V1.1

## Elaborative Design, Enhanced Usability

Newest IML craftwork makes no crack between screen panel and shell to avoid liquid ingress.

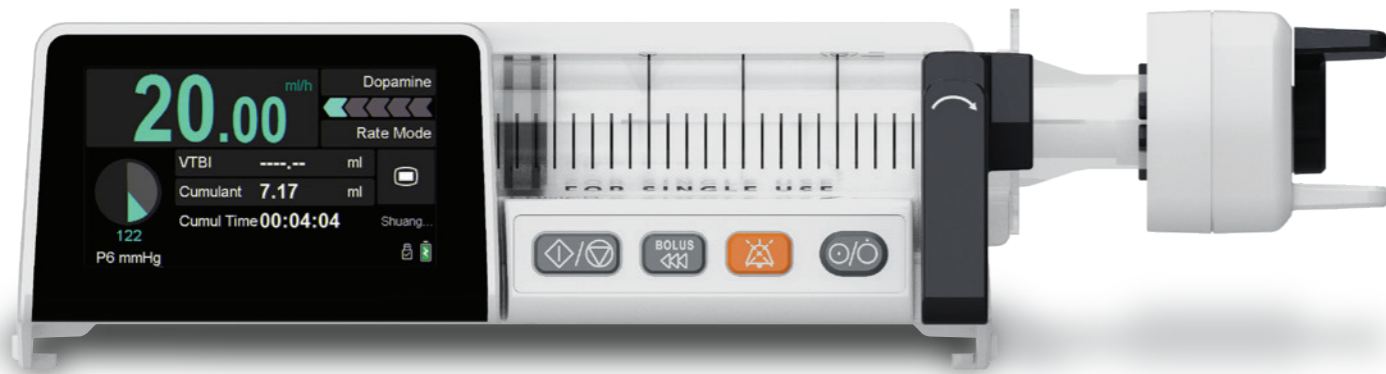

## Precise Injection



High performance motor and CPUs, ensure long-term infusion accuracy and flow rate stability.



Accuracy calibration  $\pm 1.8\%$


## Multiple Injection Modes

- The M800 was designed to perfectly meet the needs of patients requiring micro injections for a wide array of drugs. M800 would lighten caregivers' workload with up to 9 unique injection modes to choose from.

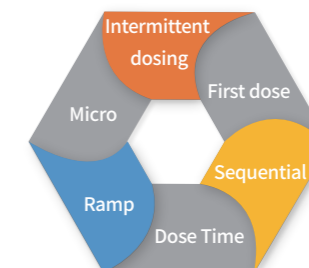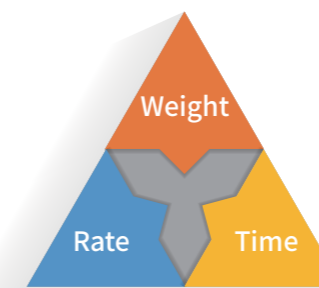


### Exceptional Craftwork

Seamless matching among parts, higher liquid ingress protection

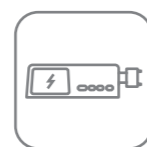

### Auto-dimming system

Ensures the proper screen brightness all the time

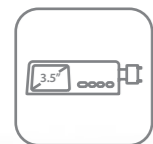

### 3.5" touchscreen

Brings effortless experience to user



### Auto-identification of syringes

Supporting 2ml, 3ml, 5ml, 10ml, 20ml, 30ml, 50ml/60ml syringe.

## Stronger combination, Flexible transport

Smaller, lighter and special handle design for various transport demands.

The modular, plug-and-play design brings the greatest flexibility of combination for various infusion scenarios.



## Intelli.™ Plan

- Comen M800 comes bundled with a drug library that can store up to 5,000 drugs to allow caregivers the flexibility of customizing drug parameters for use in different treatment plans.
- Frequently used treatment plans can be saved for future use in a variety of clinical situations. With just a few taps, the M800 will automatically load any previously saved infusion parameters based on treatment plan, increasing work efficiency while bringing to caregivers the convenience and ease they've always needed.

Add	Delete	Treatment Plan
		Treatment Plan 1--Atorvastatin--Rate Mode
		Treatment Plan 2--Adrenaline--Time Mode
		Treatment Plan 3--Atorvastatin--Rate Mode
		Treatment Plan 4--Verapamil--Time Mode

Drug	Atorvastatin	Rate Mode
Rate	200.00	ml/h
VTBI	100.00	ml
Time	00:30:00	h:m:s
	Save Treat.	